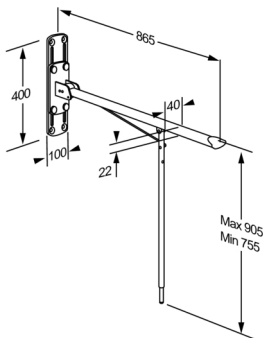




PRESSALIT®

Product Item number	VALUE support arm, height adjustable R1111000
Colour	VVS:615722410 NRF:6191644 RSK:7972413 EAN:5708590176437 HMI: 25414 NBS:N13:436 FDA:ILS/D135866 White
Product description	<div>  </div> <p>Support arm with leg, height adjustable 150 mm. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p> <p>Support arm with leg, height adjustable 150 mm. Incl. fixing bolts for brick wall. Plugs depend on wall type and must be ordered separately.</p>
Adjustability	Height adjustable
Adjustment height	150 mm
Product capacity	Maximum load: 135 kg. When the support arm is loaded with 135 kg, each wall fastener is loaded with 29 kg.
Materials	The support arm is manufactured from steel with a Rilsan coating, which is comfortable to the touch, easy to clean, and protects the metal against corrosion. The end caps are manufactured from nylon (PA).
Ergonomy	Provides firm support when you are sitting, standing, or transferring to and from a toilet. Good, broad support surface for the hand and arm. The rounded form gives a firm and comfortable grip.



PRESSALIT®

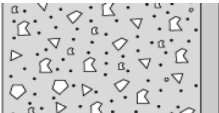
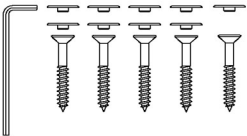

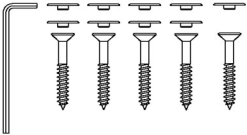
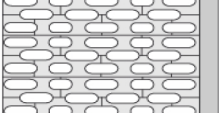
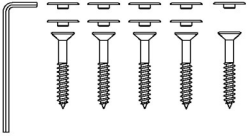
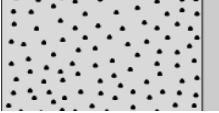
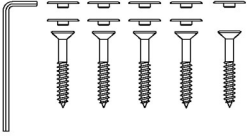
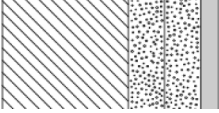
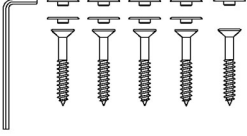
Product	VALUE support arm, height adjustable
Item number	R1111000
Features	<p>The support arm can fold up when not in use or to make more space.</p> <p>With leg for extra support.</p> <p>Height adjustable 150 mm.</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>VALUE support arms.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>For surface disinfection, chlorine-based disinfectants may be used according to the product specification and safety precautions for these products.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Klorilli</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
Declaration of conformity	CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices.



PRESSALIT®

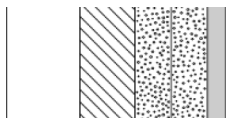
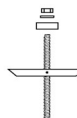


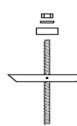

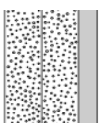

Product	VALUE support arm, height adjustable
Item number	R1111000

Installation in different wall types

	<p>Concrete >B25</p> <p>V8592 - 5 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> 
	<p>Solid Brick >= Mz12</p> <p>V8592 - 5 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> 
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8592 - 5 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> 
	<p>Aerated Concrete P4</p> <p>V8592 - 5 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> 
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8592 - 5 countersunk screws Ø10 x 80, with white screw covers are included. Plugs to be ordered separately.</p> 



PRESSALIT®

Product Item number	VALUE support arm, height adjustable R1111000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (5 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com

Pressalit A/S
 Pressalitvej 1
 DK-8680 Ry
 Tel.: +45 8788 8788
 Fax: +45 8788 8789
 E-mail: pressalit@pressalit.com

PRESSALIT®