




PRESSALIT®

Product Item number	CareSystem shower chair 410 with aperture, manually height and sideways adjustable R7568000
Colour Product description	<div>  </div> White <p>Shower chair 410 with aperture, "soft feel", with back rest and arm supports, seat height adjustable 330 mm, for mounting on horizontal track. Incl. mounting kit for brick wall. Discontinued as per 31/12-2015</p> <p>Shower chair 410 with aperture, "soft feel", with back rest and arm supports, seat height adjustable 330 mm, for mounting on horizontal track. Incl. mounting kit for brick wall. Discontinued as per 31/12-2015</p>
Adjustability Product capacity	Height and sideways adjustable Max. vertical load: 150 kg.
Materials	<p>Basic material: 20 micrometre anodized aluminium</p> <p>Seat/backrest: Polyurethane foam with marine ply core</p> <p>Seat/backrest surface: Gum lac</p> <p>Armrest surface: Two-component plastic (PVC)</p> <p>End cap, armrest: Arnite</p> <p>Lever: Stainless steel w/nylon knob</p> <p>Front plate/end covers, track: Polystyrene plastic (PS)</p>

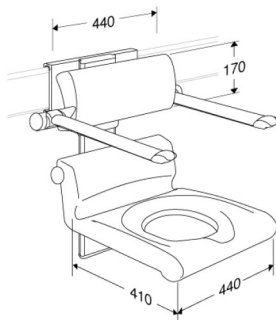

PRESSALIT®

Product	CareSystem shower chair 410 with aperture, manually height and sideways adjustable
Item number	R7568000
Ergonomy	<p>The armrests' oval shape gives an easy and comfortable surface to lean and push against.</p> <p>Soft and slip-resistant material that has a warm feel.</p> <p>Ergonomically designed seat.</p> <p>Ergonomically designed seat with aperture that makes intimate hygiene easier.</p>
Features	<p>The seat and armrests can fold up when not in use or to make more space.</p> <p>Height adjustable 330 mm - the height of the seat and the backrest can be adjusted separately.</p> <p>Can be pushed along the track.</p>
Cleaning	<p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com