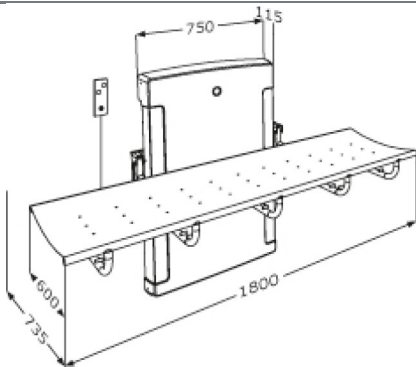




Product Item number	Nursing bench 1000, electrically height adjustable R8401000
Colour	VVS:777544180 NRF:7501953 EAN:5708590293257 HMI: 88219 NBS:N13:40,375 White
Product description	 <p>Discontinued per 31 Dec 2018 and replaced by R8403 Nursing bench 1000, powered, height adjustable 400 mm; from 480 mm - 880 mm, with wired hand control Foldable. Outer measurements 735 mm x 1800 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p> <p>Discontinued per 31 Dec 2018 and replaced by R8403 Nursing bench 1000, powered, height adjustable 400 mm; from 480 mm - 880 mm, with wired hand control Foldable. Outer measurements 735 mm x 1800 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p>
Adjustability	Height adjustable
Adjustment height	400 mm
Product capacity	Maximum load: 150 kg. When the bench is loaded with 150 kg, each mounting bolt supports 20 kg <i>This product has been tested statically with a load of 375 kg as the IEC 60601-standard requires testing with safety factor 2.5.</i>
Speed upwards	Speed up mm: 15 mm/s Speed down mm: 20 mm/s
Noise level	48 dB



Product	Nursing bench 1000, electrically height adjustable
Item number	R8401000
Materials	The bench is made from practical and durable PVC canvas which is phthalate free. The frame is powder-coated steel. Material: Steel, PVC.
Ergonomy	Pressalit Care offers a range of high quality, fixed or power adjusted nursing benches designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes. Carers can easily set the height of the bench, enabling them to obtain the correct working position. This ensures safe showering and changing for the user as well as minimizing the risk of back injury for the carer.
Features	Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall. The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space in front of the unit. Splash curtain (optional) for the nursing bench protects the carer against splashing.
Mains connection (V)	230
Actuator (motor)	LA31 Careline, (Linak), 4000 N, IPX6
Control box	CB9 Careline, (Linak), IPX6
Wired remote control	HB3 (Linak), IP66
Power consumption (maximum load)	Power consumption by maximum load, watt: 70 W Power consumption by standby: 2 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



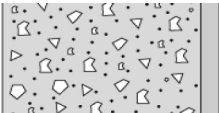
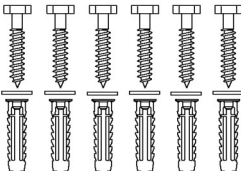

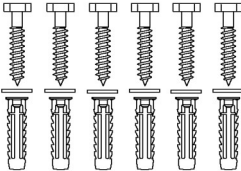
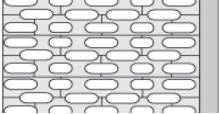
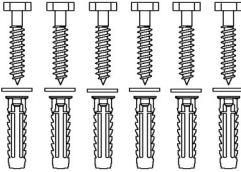
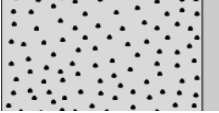
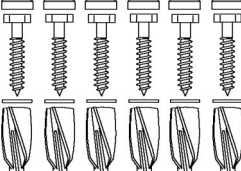
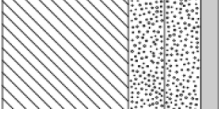
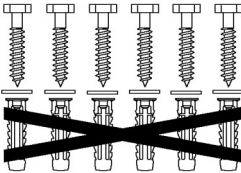
PRESSALIT®

Product	Nursing bench 1000, electrically height adjustable
Item number	R8401000
Cleaning	<p>SCT 1000 & 2000.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Accel PREvention Wipes Lingettes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>

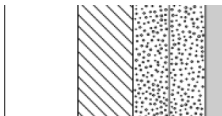
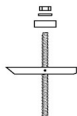

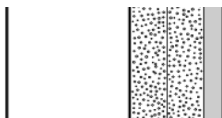
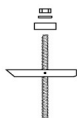

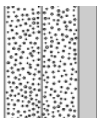



Product	Nursing bench 1000, electrically height adjustable
Item number	R8401000

Installation in different wall types

	<p>Concrete >B25</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	



Product Item number	Nursing bench 1000, electrically height adjustable R8401000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.