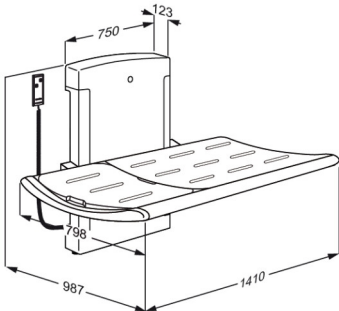




Product Item number	Nursing bench 3000, electrically height adjustable R8513021
Colour	VVS:777540234 NRF:7501901 EAN:5708590265575 HMI: 42881 NBS:N13:40,375 Lead grey
Product description	 <p>Discontinued per 31 Dec 2018 and replaced by R8533 Nursing bench 3000, powered, height adjustable 700 mm; from 300 mm - 1000 mm, with wired hand control Foldable. Incl. water collection tray with integrated outlet. Surface length 1300 mm, outer measurements 990 mm x 1410 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p> <p>Discontinued per 31 Dec 2018 and replaced by R8533 Nursing bench 3000, powered, height adjustable 700 mm; from 300 mm - 1000 mm, with wired hand control Foldable. Incl. water collection tray with integrated outlet. Surface length 1300 mm, outer measurements 990 mm x 1410 mm. Excl. mounting kit. See Technical datasheet (pdf) for information on fixings suitable for various wall types.</p>
Adjustability	Height adjustable
Adjustment height	700 mm
Product capacity	Maximum load: 200 kg. When the bench is loaded with 200 kg, each mounting bolt supports 26 kg <i>This product has been tested statically with a load of 500 kg as the IEC 60601-standard requires testing with safety factor 2.5.</i>
Speed upwards	Speed up mm: 15 mm/s Speed down mm: 20 mm/s
Noise level	48 dB



Product	Nursing bench 3000, electrically height adjustable
Item number	R8513021
Materials	Aluminium, polyurethane foam, plastics, powder-laquered steel
Ergonomy	<p>Pressalit Care offers a range of high quality, fixed or power adjusted nursing benches designed to provide a stable platform suitable for changing adults and children in a variety of settings: schools, day care centres, institutions, hospitals and private homes.</p> <p>The nursing bench 3000 series provides for maximum safety and comfort. Besides an award winning design, the series includes comfortable and durable surfaces, including polyurethane foam and aluminum.</p> <p>The flexibility of the nursing bench provides for optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p>
Features	<p>A water collection tray with integrated water outlet prevents water from splashing on carers or the floor. The flexible hose attached to the centre mounted outlet of the water collection tray, can be connected to either an outlet in the floor or on the wall. Adjustable neck support which can be moved to either end of the 1800 mm nursing benches.</p> <p>The slightly curved form makes the bench extremely comfortable and safe.</p> <p>The integrated gas cylinders enable the bench to be folded out of the way when not in use to maximise floor space.</p>
Mains connection (V)	230
Actuator (motor)	LA31 Careline, (Linak), 6000 N, IPX6
Control box	CB9 Careline, (Linak), IPX6
Wired remote control	HB3 (Linak), IP66
Power consumption (maximum load)	<p>Power consumption by maximum load, watt: 70 W</p> <p>Power consumption by standby: 2 W</p>
Warranty	<p>Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.</p>



Product	Nursing bench 3000, electrically height adjustable
Item number	R8513021
Cleaning	<p>SCT 3000.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes Universal Wipes (Clinell) Sporidical Wipes (Clinell) CHLOR-CLEAN Tablets (chemform) WipeClean Ethanol Disinfection 80% (PLUM) Wetwipe Ethanol Disinfection 70% (wetwipe)</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>



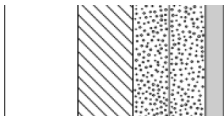
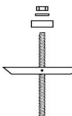

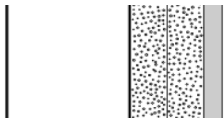
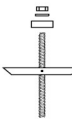

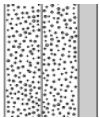

PRESSALIT®

Product	Nursing bench 3000, electrically height adjustable
Item number	R8513021

Installation in different wall types

	<p>Concrete >B25</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. To be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8658 - Mounting kit, 6 coach screws Ø10 x 70 - plugs GB14. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8608 - Mounting kit, 6 coach screws Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.</p>	



Product Item number	Nursing bench 3000, electrically height adjustable R8513021		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.