




<b>Product</b> <b>Item number</b>	PLUS wash basin bracket with cordless remote control, electrically height adjustable and sideways adjustable R4782
<b>Product description</b>	<p>           VVS: 654003230            NRF: 6132539            RSK: 7972597            EAN: 5708590309354            HMI: 71754            NBS: N13:30,335         </p>  <p>           PLUS wash basin brackets are available in electrically, pneumatically and manually height adjustable versions. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.         </p> <p>           PLUS wash basin bracket, powered with cordless remote control, electrically height adjustable by 300 mm and sideways adjustable. Recommendation: Flexible feed and waste system R2064/R2065. (UK markets R2074/R2075).            Materials: Aluminium and ABS polymer.            Incl. fixing bolts for wash basin. For fitting to horizontal PLUS wall track (not included).            Will be discontinued no later than 12/31-2019. Subject to availability.         </p>
<b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b>	<p>Height and sideways adjustable</p> <p>300 mm</p> <p>           Maximum load: 100 kg.            Suitable for wash basins that weigh 16-22 kg.            When the wash basin is loaded with 100 kg, each wall fastener will be carrying 22 kg.  <i>Wash basin brackets need to be able to carry the weight of the wash basin, as well as support the weight of the person who is leaning against the wash basin or who is using the wash basin to rise from seated position.</i>  <b>Neither wash basin nor wash basin brackets are intended to carry a person, and the test methods applicable to wash basin brackets reflect this: ISO 17966 prescribes a static test on the front edge of the wash basin at 75% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 100 kg and the wash basin weighs 20 kg, it is – according to ISO 17966 – the equivalent of a user weighing 106 kg leaning on the front edge of the wash basin.</b> </p>
<b>Speed upwards</b>	<p>           Speed up mm: 20 mm/s            Speed down mm: 23 mm/s         </p>

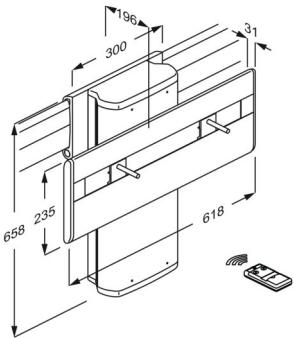

# PRESSALIT®

<b>Product</b>	PLUS wash basin bracket with cordless remote control, electrically height adjustable and sideways adjustable
<b>Item number</b>	R4782
<b>Noise level</b>	45 dB
<b>Materials</b>	The wash basin bracket is primarily made of 10-micrometre anodized aluminium. The ends of the bracket cover are capped with ABS polymer. The bracket is available with a cordless remote control made of ABS polymer.
<b>Ergonomy</b>	The height of the wash basin is easily adjusted by means of the cordless remote control. The wash basin bracket suits all types of bolt-mounted wash basin. Its rounded, sleek design greatly facilitates cleaning.
<b>Features</b>	An independant user will be able to make the necessary adjustments to an electrically operated wash basin bracket themselves without any problems.
<b>Mains connection (V)</b>	230
<b>Actuator (motor)</b>	LA27CS (Linak) 1500N, IPX6
<b>Control box</b>	TR6 (Linak) IPX0
<b>Wired remote control</b>	SEAV (Linak) IP54
<b>Power consumption (maximum load)</b>	Power consumption by maximum load, watt: 75 W Power consumption by standby: 1,5 W
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS wash basin bracket with cordless remote control, electrically height adjustable and sideways adjustable  R4782
<b>Cleaning</b>	<p>PLUS &amp; SELECT wash basin brackets. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Clinimax Difficil-S Erisan OXY+</p> <p><b>NON-DISINFECTANTS:</b> Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

[pressalit.com](http://pressalit.com)

## PRESSALIT®