




PRESSALIT®

Product	PLUS backrest, height and sideways adjustable
Item number	R6530112000
Colour	<p>VVS:777538600 NRF:6190622 EAN:5708590314488 HMI: 71808 NBS:N13:426 FDA:INN/D135865 White</p> 
Product description	<p>PLUS backrest 385 mm, soft feel. Manually height adjustable by 150 mm and sideways adjustable. Anthracite cover plate for vertical wall track. Materials: Aluminium, polyurethane foam with a plywood core. For fitting to horizontal PLUS wall track (not included).</p> <p>PLUS backrest 385 mm, soft feel. Manually height adjustable by 150 mm and sideways adjustable. Anthracite cover plate for vertical wall track. Materials: Aluminium, polyurethane foam with a plywood core. For fitting to horizontal PLUS wall track (not included).</p>
Adjustability	Height and sideways adjustable
Adjustment height	150 mm
Product capacity	100 kg.
Materials	The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer.
Ergonomy	Soft and slip-resistant material that has a warm feel.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.

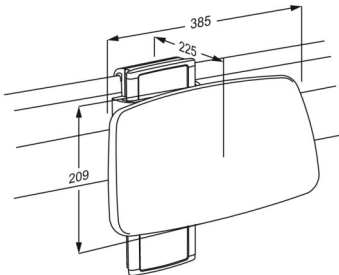

PRESSALIT®

Product	PLUS backrest, height and sideways adjustable
Item number	R6530112000
Cleaning	<p>PLUS backrest.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®