




# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS shower seat 450 with aperture, manually height and sideways adjustable R7451112000
<b>Colour</b>	VVS: 777536800 NRF: 6311194 EAN: 5708590312774 HMI: 71795 NBS: N13:448 FDA: INN/D135865 White
<b>Product description</b>	 <p>PLUS shower seat 450 with aperture, with back- and armrests. Manually height adjustable by 175 mm and sideways adjustable. Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track.</p> <p>Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS. For fitting to horizontal PLUS wall track (not included).</p> <p>PLUS shower seat 450 with aperture, with back- and armrests. Manually height adjustable by 175 mm and sideways adjustable. Armrests are independently height adjustable. The seat and the armrests are foldable. Anthracite cover plate for vertical wall track.</p> <p>Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. The non-slip gripping surface of the armrests are made of soft SEBS. For fitting to horizontal PLUS wall track (not included).</p>
<b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b>	Height and sideways adjustable 175 mm Maximum user mass (seat): 200 kg. Maximum load for armrests is 55 kg. When the seat is loaded with 200 kg, each wall fastener is loaded with 23 kg. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>

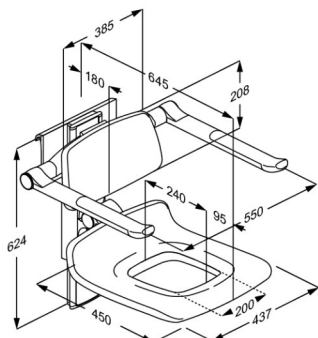

# PRESSALIT®

<b>Product</b>	PLUS shower seat 450 with aperture, manually height and sideways adjustable
<b>Item number</b>	R7451112000
<b>Materials</b>	The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. The arms are primarily constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.
<b>Ergonomy</b>	Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel. The oval shape of the armrests gives an easy and comfortable surface to lean on and push against.
<b>Features</b>	The seat, backrest and armrests can be adjusted simultaneously thanks to a safe and easy click system. The armrests can also be adjusted independently of the seat and backrest. The seat and armrests can fold up when not in use or to maximise floor space.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



<b>Product</b>	PLUS shower seat 450 with aperture, manually height and sideways adjustable
<b>Item number</b>	R7451112000
<b>Cleaning</b>	<p>PLUS shower seat.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b> Erisan OXY+ Tvättspit 70% CaviCide surface disinfectant decontaminant cleaner Rodalon Desinfektionsspray Terralin protect (SCHÜLKE+) Dismonzon Plus (BODE Chemie GmbH)</p> <p><b>NON-DISINFECTANTS:</b> Wetrok Caledor (Wetrok) Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>



# PRESSALIT®

<b>Product</b>	PLUS shower seat 450 with aperture, manually height and sideways adjustable
<b>Item number</b>	R7451112000

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®