




Product Item number	SELECT toilet lifter, manually height adjustable with crank handle, for wall outlet R8125000
Colour Product description	<div> <div>White</div>  </div> <p>Discontinued per 15 May 2017 and replaced by R8122. Wallhung toilet 540 mm to be ordered separately (R2091).</p> <p>Toilet lifter, 300 mm height adjustable, with crank handle. Incl. cistern with dual flushing system preset at 3/6 litres, wall-hung toilet pan (pan projection 540 mm), and R37 white Dania toilet seat with extra strong cross-bar hinge D92. The integral vertical side profiles allow for addition of PLUS support arms. For wall-outlet. Materials: Anodised aluminium and polystyrene (PS). Incl. mounting kit for brick wall.</p> <p>Discontinued per 15 May 2017 and replaced by R8122. Wallhung toilet 540 mm to be ordered separately (R2091).</p> <p>Toilet lifter, 300 mm height adjustable, with crank handle. Incl. cistern with dual flushing system preset at 3/6 litres, wall-hung toilet pan (pan projection 540 mm), and R37 white Dania toilet seat with extra strong cross-bar hinge D92. The integral vertical side profiles allow for addition of PLUS support arms. For wall-outlet. Materials: Anodised aluminium and polystyrene (PS). Incl. mounting kit for brick wall.</p>
Adjustability Adjustment height Product capacity	Height adjustable 300 mm Maximum load: 200 kg
Materials Ergonomy	<p>The toilet lifter is mainly constructed of anodised aluminium. The front plate and the bottom plate are made of polystyrene (PS).</p> <p>Toilet lifters are available in three different colours which all harmonise with modern bathrooms, resulting in a comfortable and inviting toilet environment. A minimum of grooves, slots and holes makes cleaning fast and effective, and also results in an attractive and exclusive design.</p>

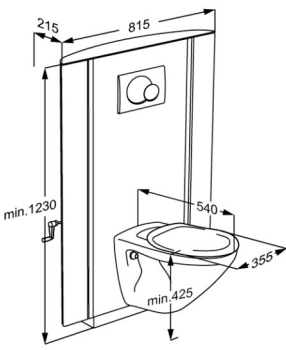

PRESSALIT®

Product	SELECT toilet lifter, manually height adjustable with crank handle, for wall outlet
Item number	R8125000
Features	<p>The toilet can be raised and lowered using the crank handle. The crank handle can be positioned on either side of the lifter – without the need for tools.</p> <p>It is possible to adjust the height of the toilet lifter using an electric drill as an alternative to the crank handle.</p> <p>The toilet can be adjusted 300 mm, from 425 mm to 725 mm in height , so the height can always be perfectly matched to mobility aids such as a wheelchair.</p>
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.

pressalit.com