




Product Item number	SELECT toilet lifter, manually height adjustable with crank handle, for floor outlet R8117000
Colour Product description	<div>  </div> <p>White</p> <p>Discontinued per 9 Jan 2017 and replaced by R8112. Wallhung toilet 700 mm to be ordered separately (R2099).</p> <p>Toilet lifter, 300 mm height adjustable, with crank handle. Incl. cistern with dual flushing system preset at 3/6 litres, wall-hung toilet pan (pan projection 700 mm), and R27 white Dania toilet seat with extra strong cross-bar hinge D92. The integral vertical side profiles allow for addition of PLUS support arms. For floor-outlet. Materials: Anodised aluminium and polystyrene (PS). Incl. mounting kit for brick wall.</p> <p>Discontinued per 9 Jan 2017 and replaced by R8112. Wallhung toilet 700 mm to be ordered separately (R2099).</p> <p>Toilet lifter, 300 mm height adjustable, with crank handle. Incl. cistern with dual flushing system preset at 3/6 litres, wall-hung toilet pan (pan projection 700 mm), and R27 white Dania toilet seat with extra strong cross-bar hinge D92. The integral vertical side profiles allow for addition of PLUS support arms. For floor-outlet. Materials: Anodised aluminium and polystyrene (PS). Incl. mounting kit for brick wall.</p>
Adjustability Adjustment height Product capacity	Height adjustable 300 mm Maximum load: 200 kg
Materials Ergonomy	<p>The toilet lifter is mainly constructed of anodised aluminium. The front plate and the bottom plate are made of polystyrene (PS).</p> <p>Toilet lifters are available in three different colours which all harmonise with modern bathrooms, resulting in a comfortable and inviting toilet environment. A minimum of grooves, slots and holes makes cleaning fast and effective, and also results in an attractive and exclusive design.</p>

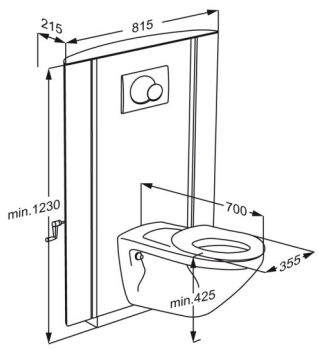

PRESSALIT®

Product	SELECT toilet lifter, manually height adjustable with crank handle, for floor outlet
Item number	R8117000
Features	<p>The toilet can be raised and lowered using the crank handle. The crank handle can be positioned on either side of the lifter – without the need for tools.</p> <p>It is possible to adjust the height of the toilet lifter using an electric drill as an alternative to the crank handle.</p> <p>The toilet can be adjusted 300 mm, from 425 mm to 725 mm in height , so the height can always be perfectly matched to mobility aids such as a wheelchair.</p>
Cleaning	<p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®