




Product Item number	PLUS replacement shower seat 450, manually height adjustable R7402000
Colour Product description	White  PLUS shower seat 450, with backrest. Manually height adjustable by 360 mm. The seat and the armrests are foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For fitting to vertical PLUS wall track (not included). PLUS shower seat 450, with backrest. Manually height adjustable by 360 mm. The seat and the armrests are foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For fitting to vertical PLUS wall track (not included).
Adjustability Adjustment height Product capacity	Height adjustable 360 mm Maximum user mass: 200 kg. When the seat is loaded with 200 kg, each wall fastener is loaded with 24 kg. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
Materials Ergonomy Features	The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS. Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel. The seat can fold up when not in use or to maximise floor space. The seat and backrest can be adjusted simultaneously thanks to a safe and easy click system.

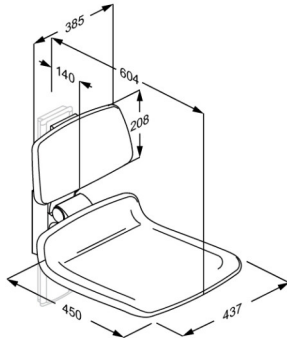

PRESSALIT®

Product	PLUS replacement shower seat 450, manually height adjustable
Item number	R7402000
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
Cleaning	<p>PLUS shower seat.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®